i

See e

K

breo

breo

CONTENTS

Safety Information **03**

Supplied Items **08**

Features **09**

Instructions for Use **10**

Uninstalling Batteries **18**

Trouble Shooting **20**

Care and Storage **21**

Guarantee/Warranty **24**

Product Specification **26**

EN**Precautions**

Please read this instruction manual carefully before using. The company assumes no responsibility for any consequences caused by misuse or failure of following these instructions.

	Do not use this product while driving.
	Please remove any accessories you might be wearing or other objects that might affect the massaging process before use.
	Do not use the product in the bathroom or other humid conditions.
	Keep the product away from water, fire, and corrosive environment.
	Please comply with relevant safety regulations while using the product under specific circumstances.
	Do not use the product on surgical areas. Do not use the product on areas of bruises, damaged skin, swelling, or inflammation.

	<p>Do not use the product if you have cardiovascular (or once had), acute, or dermatologic diseases, malignant tumors, or other abnormal symptoms.</p> <p>Do not use the product in case of hematologic disorders with hemorrhagic symptoms, osteoporosis, soft tissue diseases, or strokes.</p> <p>In case of restricted physical capabilities, physical or mental problems, or mental disabilities, do not use the device without medical consultation first.</p> <p>Stop the massage if you find it unpleasant.</p>
	<p>Do not use the device during menstruation or pregnancy.</p>
	<p>This device is not intended for use by children under the age of 6. Children aged between 8 and 12 must be supervised while using the device.</p>
	<p>Clean with alcohol wipes. Avoid contact with water.</p>
	<p>Avoid contact with sharp objects.</p>

EN

The device is for private use only and cannot be used as medical equipment. It is not intended for medical diagnosis, treatment, or commercial purposes.

Please use only the supplied charging device to charge this product and ensure electrical safety.

If any fault or malfunction should occur during use, cut off the power immediately. Never attempt to disassemble or repair the device yourself.

About this device

※ Usage Frequency

Use it in the morning and evening for no more than 10 minutes each time.

※ Before you use this device,

Please note

- 1 Do not use the device with wet hands.
- 2 Remove your contact lenses, hair clips, hair accessories, etc. to avoid discomfort during the massage.
- 3 If you feel pain or burning in or around your eyes, please discontinue using the device and consult a physician.

※ How to use this device

When using this device, please keep your eyes closed, lean back in the chair and completely relax your body.

Remove your contact lenses, glasses, false eyelashes, etc.


EN ※ Occurrences after using this device

Redness may occur around your eyes, but they will return to normal after a short time.

Indentations will appear on your face, but they will disappear after a short time.

※Additional info

The "buzzing", "popping" and "booming" sounds heard during the massage are produced by fine-tuning rather than the malfunction of the product.

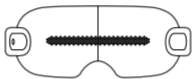iSee K

charging cable

instruction manual

EN

PRODUCT FEATURES


**Vibrational
Massage**


**Meditation
Guidance**

Visualization

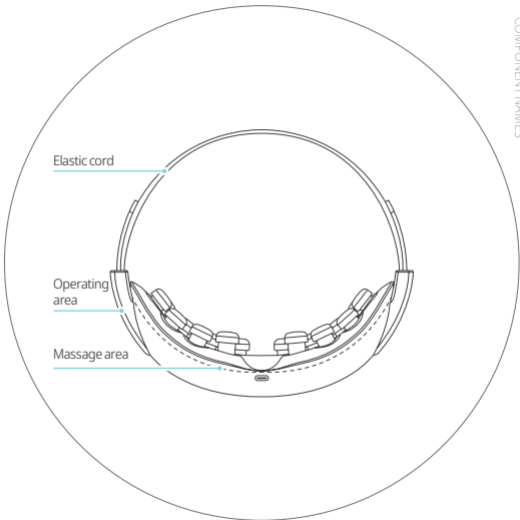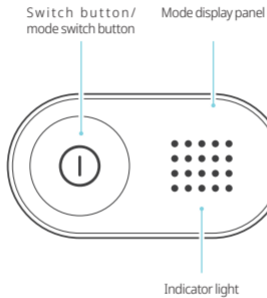Charging

1 Connect the USB to the port on the massager.

2 Put the adapter in the socket.

After a "beep", the charge indicator flashes on the control panel.

The indicator lights up continuously when the device is fully charged.

※ The massager and the adapter may become warm during normal charging process.

Time for charging

1 Charging takes up to 2-3 hours.

2 Once fully charged, the massager can be used for 10 times (10 minutes a time).


(Use time might differ when operating in different temperatures.)

3 When the indicator flashes, charge immediately.

Fully Charged

Unplug the adapter and remove the USB from the massager.

EN

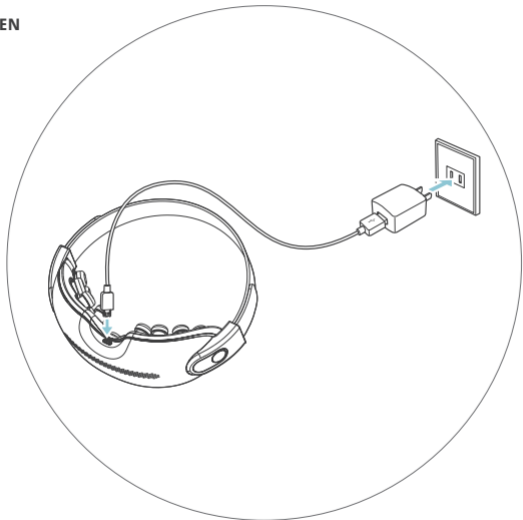How to wear this device

- 1 Pull out the elastic cord of this device.
- 2 Put the device around your head and adjust it to the appropriate angle.
- 3 Press and hold the switch button for 2 seconds to start the device. The massage will begin after the voice prompt is heard.

※ Please be careful not to wear it too tightly.

Please make sure that you have removed your contact lenses, glasses, false eyelashes, etc. before using this device.

EN

Correct usage

The device cannot be used while being charged. Before you use the device, please unplug the power adapter after charging is completed.

1 Press and hold the switch button for 2 seconds to start up. The massage procedure will begin after the voice prompt. The eye exercises mode is enabled by default with the indicator light staying on.

2 Tap the switch button to switch modes if desired.

Mode 1:Eye exercises mode

Mode 2:Beauty mode

Mode 3:Eye protection mode

Mode 4:Touching mode

Mode 5:Music-driven mode

3 Music playback

Turn on the Bluetooth of your cell phone, and search "iSee K"; After successful connection, the massager is ready to play music of your cell phone; To customize operating pattern of the nodes with your music, switch to Model 5 and raise the volume of your cell phone.

4 Description of Power Indicator

State	Power indicator display
Normal Operation (Bluetooth Disconnected)	Flashing Purple Light
Normal Operation (Bluetooth connected)	Solid Purple Light
Low Battery	Flashing Red Light
Charging	Flashing Red Light
Fully Charged	Solid Purple Light
Overvoltage Charging	No light

5 Wear this device around your head. When you are using this device, please keep your eyes closed, lean back in the chair and completely relax your body.

6 A voice prompt will be heard after operating for 10 minutes. Then the indicator light will go off and the device will automatically shut down.

7 How to cut off power manually

When you press and hold the switch button for over 2

seconds, the device will automatically shut down.

8 Please remove the device from around your head after you finish using it.

※ Your vision will be blurred for a short time after the eye massage. It will return to normal after you keep them closed for a short while.

Operation steps for battery removal

- 1 Cut off the external power supply, and remove the product housing;
- 2 Disconnect the power supply to the battery;
- 3 Remove the battery;
- 4 Place the waste battery at a fixed recycling station;
- 5 Please observe local environmental protection laws and waste disposal guidance and protect the environment.

EN

Notes

- It must only be supplied at safety extra low voltage corresponding to the marking on the appliance.
- This appliance contains batteries that are non-replaceable.
- Remove the battery out of the appliance before scrapping the appliance;
- The appliance must be disconnected from power supply before removing the battery;
- The battery should be disposed of safely.
- This appliance can be used by children aged from 8 years and above and persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instruction concerning use of the appliance in a safe way and understand the hazards involved. Children shall not play with the appliance. Cleaning and user maintenance shall not be made by children without supervision.

The device may only be repaired by qualified personnel!

Troubleshooting

Please perform the following checks in case of any failure occurring during operation. If the failure cannot be eliminated or the product still cannot work ideally, please cut off power supply and ask a professional technician to repair it. Never attempt to disassemble the product.

Symptom	Possible cause	Basic solution
Power-on failure	Check whether the battery level is too low	If yes, please charge the battery timely
Low massage intensity	<ul style="list-style-type: none"> ① check whether the batter level is low ② Please check whether the volumn of your device gets to low 	<ul style="list-style-type: none"> ① If yes, please charge the battery timely ② Please raise the volumn if so.

If the problem still cannot be solved after troubleshooting, please contact us at service_global@breo.com

EN

Storage

When the device has not been used for over 3 months, take out and recharge it before use. Do not place it next to a wall or on the ground or rooftop.

Make sure the room is well-ventilated. Keep the massager away from direct sunlight and any corrosive gas.

Travel

The shockproof packaging can be used for shipping by air, rail, road, and sea. The device should be protected from rainy and snowy weathers. Avoid placing it upside down and crashing it.

Power Check

Check the voltage rating before switching on the device to ensure the required voltage match the available power source.

Maintenance

Keep the device dry, clean, and away from rats. To clean the device, please make sure the power is off. Put the device in a safe place after use. Do not place anything heavy on it. Always bring your product to a qualified professional for any service or repair. Do not disassemble the product yourself.

Disposal and Recycling

Do not throw away the product with the normal household waste at the end of its life, but follow the local environmental regulations and guidance for disposal.

FCC Regulatory Compliance

This device complies with Part 15 of the FCC Rules. Operation is subject to the following two conditions: (1) this device may not cause harmful interference, and (2) this device must accept any interference received, including interference that may cause undesired operation.

The statements should be displayed in the user manual: changes or modifications not expressly approved by the party responsible for compliance could void the user's authority to operate the equipment.

This equipment has been tested and found to comply with the limits for a Class B digital device, pursuant to Part 15 of the FCC Rules. These limits are designed to provide reasonable protection against harmful interference in a residential installation. This equipment generates, uses and

EN

can radiate radio frequency energy and, if not installed and used in accordance with the instructions, may cause harmful interference to radio communications. However, there is no guarantee that interference will not occur in a particular installation.

If this equipment does cause harmful interference to radio or television reception, which can be determined by turning the equipment off and on, the user is encouraged to try to correct the interference by one or more of the following measures:

- Reorient or relocate the receiving antenna.
- Increase the separation between the equipment and receiver.
- Connect the equipment into an outlet on a circuit different from that to which the receiver is connected.
- Consult the dealer or an experienced radio/TV technician for help.

This equipment complies with FCC radiation exposure limits set forth for an uncontrolled environment. This transmitter must not be co-located or operating in conjunction with any other antenna or transmitter.

Guarantee/Warranty:

The warranty period is 2 years from the date of purchase. This includes a 12-month guarantee from the date of purchase. This guarantee is valid for the whole of Europe. The purchase receipt is regarded as proof of purchase. The customer undertakes to read the instructions carefully before making a claim under the guarantee/warranty in order to rule out any operating errors.

The product has been thoroughly checked for quality and functionality before dispatch. Exclusion of the guarantee occurs in the following cases:

- in the event of changes to the product without our approval
- in the event of use other than for the intended purpose
- in the event of defects caused by improper maintenance and/or repair
- in the event of commercial use of the article
- in the event of damage caused by falling

When making a claim under the guarantee/warranty, please return the complete product including accessories in a packaging safe for transport to the point of sale from which the product was purchased. Describe the existing defect in the product and enclose the original proof of purchase with the shipment.

The seller reserves the right to subsequent performance before refunding the purchase price. Necessary transport costs will be reimbursed to the customer if a justified claim of the legal warranty/guarantee is established. The exercise of the statutory rights in respect of defects shall be free of charge. These rights are not limited by the guarantee.

SHP International Trading GmbH (Rathausstraße 5, 96342 Stockheim, Germany) hereby declares that this product is in conformity with the essential requirements of the applicable European Directives.

Intended use

This device (iSeek) is intended for use as a massager only. Any other use shall be deemed inappropriate.

Product name	Eye Massager
Product model	iSee K
Rated voltage	5V 
Rated power	5W
Battery capacity	1050mAh
Product weight	230g
	L196×W78×H82(mm)

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