Personenwaage 7-in-1

Personal Scale 7-in-1

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Personenweegschaal 7-in-1

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Operating instructions

Thank you for purchasing our SEECODE Personal Scale 7-in-1, which allows you to easily and quickly measure body fat, body hydration (water percentage), muscle and bone mass, BMI and calories for up to 10 people.

To ensure that you enjoy your product for a long time and can use it to the full, please read these operating instructions carefully before using the SEECODE 7-in-1 bathroom scales to avoid damage due to incorrect operation. Please pay special attention to the safety information. If you give the SEECODE Personal Scale 7-in-1 to a third party, this instruction manual should also be handed over. Please keep the operating instructions for future reference.

The scale is only suitable for home use and not for professional use. The measured values are only to be understood as reference values. Our product is only intended for use by persons over the age of 10.

Scope of delivery

Personal scale 7-in-1

Operating instructions

Technical data

Dimensions: 31 x 31 x 2,3 cm

Weight: 2 KG

Material: ABS and glass

Batteries: 4 x AAA lithium batteries (not included)

Capacity: 6kg to 180kg

Age range: 10 to 80 years

Height range: 100cm to 240cm

Weight graduation: 0.1kg

Body fat / total body water / muscle mass division: 0.1%

Bone mass: 0.1kg

Kilocalorie indication: 1Kcal.

Automatic switch-off

Low battery indicator

Possible switching between three different units: kg/lb/st.

Features of the SEECODE personal scale

This product uses advanced bioelectrical impedance technology to measure body fat percentage and hydration by passing a safe electrical impulse through the body.

The difference in electrical conductivity of body fat and other tissues is analysed and processed by the CPU (processor) to determine the body fat, hydration and fluid percentage of the body.

All results are displayed on the LCD screen.

High-precision strain gauge sensors ensure very accurate results.


The tempered glass platform of high quality and modern design makes the scale safe and visually appealing.

The memory for the data of up to 10 people offers more convenience in daily use.

The personal data is stored, not the body values, so that this information is not accessible to third parties.

Inserting the batteries

Open the lid on the back of the scale and insert 4 AAA batteries (not included). Please pay attention to the given polarity (+ / -).

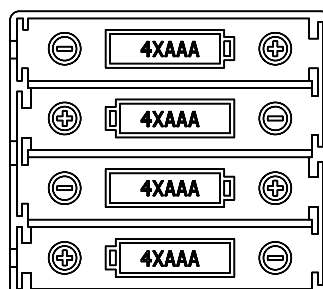Operating buttons for basic setting

 Key: Switch on, select menu and confirm

 Button: Up

 Key: Down

UNIT: Selection of units between kg/lb/st

Important instructions for operating the personal scale

Place the scale on a hard and level surface.

Remove your shoes and socks and keep the skin of your feet in contact with the 4 electrodes when measuring body fat, hydration, muscle, bone, BMI and calories.

Stand still and upright to measure your weight and body values.

The glass platform can be slippery when wet, so please ensure that both the glass and your feet are dry before use.

Your body score will be displayed on the screen within seconds. These results run consecutively and are repeated three times.

This scale is a highly accurate measuring device. Never jump or rock on the scale.

The scale switches off automatically when it is no longer in use.

Please check the battery level if the scale is not working properly. Change the batteries if necessary.

Operating the personal scales


For weight measurement only:

Switch on the scale and stand on it until the weight flashes. The scale switches off automatically when it is no longer in use.

Setting your personal data for measuring weight, body fat, hydration, muscle mass, bone mass, BMI and calories.

The personal data (man/woman, age and height) are saved once you have set them under a code number (P0 to P9). Simply select your code number when you use the scale again. Your preset information will be displayed so you do not have to program it again.

Press the „SET“ button for a few seconds to set your personal values. When the screen shows „0.0kg“ press the  or  button to assign the appropriate code number to your person. Once you have reached your code number, please press „SET“ to confirm.

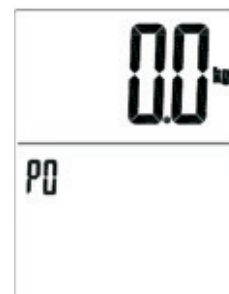Press  or  to select Man  or Woman .

When your gender symbol flashes on the screen, press „SET“ to confirm.

Press $\wedge$ or $\vee$ to set your age (10 to 80 years). When „Age“ and the number on the screen flash, press „SET“ to confirm.


Press $\wedge$ or $\vee$ to set your height (100 cm to 240 cm). When the number and „cm“ on the screen flash and your height is set in cm, then press „SET“ to confirm.

Automatic value sequence in the display

When you use the scale, all values are displayed one after the other:
Weight / Body fat / Water percentage / Muscle mass / Bone mass / BMI / KCAL.

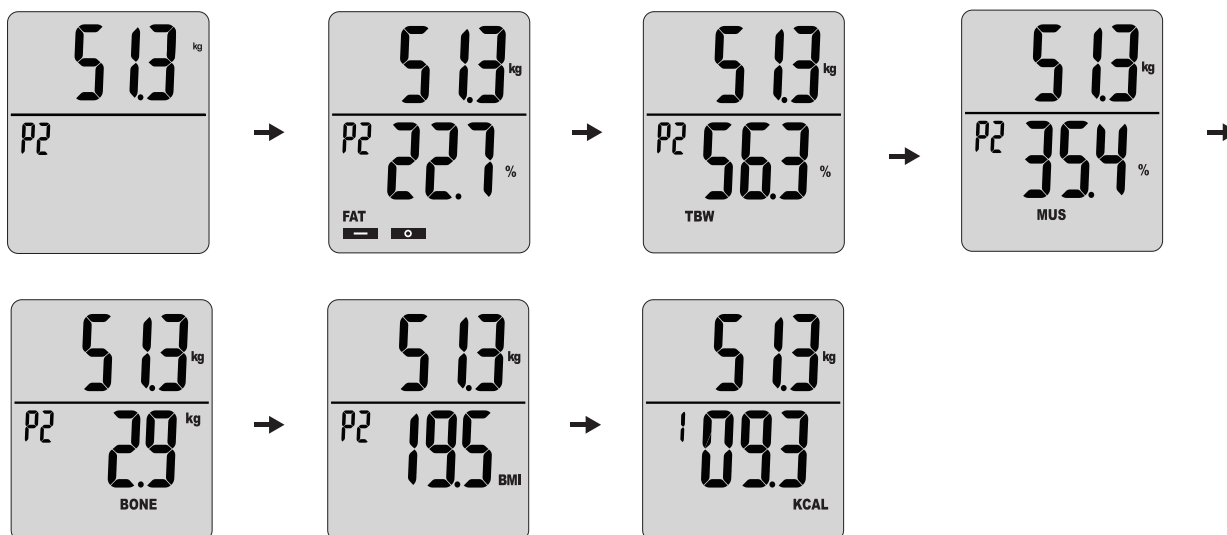Tables to classify your health values

You can verify your measured values in the following tables.

Table of body fat and hydration ranges

Age	Female		Male		Display
	Fat%	Hydration%	Fat%	Hydration%	
≤30	4.0-20.5	66.0-54.7	4.0-15.5	66.0-58.1	Low fat percentage
	20.6-25.0	54.6-51.6	15.6-20.0	58.0-55.0	Healthy fat percentage
	25.1-30.5	51.5-47.8	20.1-24.5	54.9-51.9	High fat percentage
	30.6-80.0	47.7-13.8	24.6-80.0	51.8-13.8	Overweight
≥30	4.0-25.0	66.0-51.6	4.0-19.5	66.0-55.4	Low fat percentage
	25.1-30.0	51.5-48.1	19.6-24.0	55.3-52.3	Healthy fat percentage
	30.1-35.0	48.0-44.7	24.1-28.5	52.2-49.2	High fat percentage
	35.1-80.0	44.6-13.8	28.6-80.0	49.1-13.8	Overweight

Table of muscle and bone mass

Age	Female		Male	
	Muscle mass %	Bone mass %	Muscle mass %	Bone mass %
16-50	≥31	≥8	≥38	≥10

BMI parameters:

BMI parameters			
Low	Normal	Overweight	too high fat percentage
≤18.5	18.5~25	25~30	≥ 30

Energy (Kcal):


The energy needed is consumed to maintain healthy breathing, keep warm and other basic life activities, this takes up 70% of the daily requirement.

Male			
Age	Normal weight	Basic metabolic value (Kcal/kg/day)	Recommended calorie consumption per day
15-17	59.8	27	1610
18-29	64.7	24	1550
30-49	67	22.3	1170
50-69	62.5	21.5	1110
≥70	56.7	21.5	1010

Female			
Age	Normal weight	Basic metabolic value (Kcal/kg/day)	Recommended calorie consumption per day
15-17	51.4	25.3	1300
18-29	51.2	23.6	1210
30-49	54.2	21.7	1170
50-69	53.8	20.7	1110
≥70	48.7	20.7	1010

Important information, please note

- Persons with pacemakers must not use this product.
- For persons with metal implants, the measurements may be inaccurate, the results may be falsified.
- To ensure good comparability of the values, we recommend taking the measurement at the same time every evening.
- If „Lo“ appears on the LCD display, please replace the batteries.

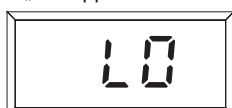The indication „ERR“ on the LCD display shows an overload of the scale.

- Further, the indication „Err“ shows that the measurement has an error.
- If the body fat percentage is less than 5% or more than 50%, the display shows „Er2“.

Care and maintenance

The surface can be cleaned with a soft, slightly damp cloth.
Do not use abrasive cleaners or harsh chemicals for cleaning.
Always keep the scale away from water and heat, e.g. sunlight.

Warranty / Guarantee

The warranty period is 24 months from the date of purchase. This includes a 6-month warranty from the date of purchase.
The territorial scope of the warranty protection is Europe-wide. The proof of purchase is valid. The product has been thoroughly checked for quality and functionality before shipment. An exclusion of the guarantee takes place under the following cases:

- in case of modifications of the system without our approval
- in case of use other than the above mentioned
- in case of damage caused by falling

In the event of a claim under the guarantee/warranty, please return the complete product in transport-safe packaging to the point of sale from which the product was purchased. Describe the detected defect in the product and include the original proof of purchase with the shipment. The seller reserves the right to remedy the defect before refunding the purchase price. Necessary transportation costs will be reimbursed to the clientele if a justified claim under the statutory warranty/guarantee is established. The use of the legal rights for defects is free of charge.

CE Declaration of Conformity

Mobiset GmbH hereby declares that this product is in conformity with the essential requirements of the applicable European Directives.

General information

Copyright

This document is protected by copyright. Any duplication or reprinting, including excerpts, as well as the reproduction of illustrations, even in modified form, is only permitted with the written consent of the manufacturer.

Notes on environmental protection

 The packaging materials used can be recycled. Dispose of packaging materials that are no longer required in accordance with the locally applicable regulations.

 This product must not be disposed of in normal household waste within the European Union. Dispose of the device in accordance with the regional regulations.

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